

Product information

Information about other products is available at: www.demeditec.com



User's Manual

Anti-LC-1 ELISA

Enzyme immunoassay for the determination of IgG antibodies to formiminotransferase-cyclodeaminase in human serum and plasma



DE4276



96

***Please use only the valid version of the package insert provided with the kit.
Verwenden Sie nur die jeweils gültige, im Testkit enthaltene, Arbeitsanleitung.
Si prega di usare la versione valida dell'inserto del pacco a disposizione con il kit.
Por favor, se usa solo la version valida de la metodico técnico incluido aqui en el kit.***

Table of Contents / Inhaltsverzeichnis / Tabella die Contenuti / Tabla de Contenidos

1	INTENDED USE	3
2	PRINCIPLE OF THE TEST	3
3	PATIENT SAMPLES.....	3
4	TEST COMPONENTS FOR 96 DETERMINATIONS.....	4
5	ASSAY PROCEDURE	5
6	DATA PROCESSING	5
7	REFERENCE VALUES	6
8	CHARACTERISTIC ASSAY DATA.....	6
9	INCUBATION SCHEME	7
10	SAFETY PRECAUTIONS.....	8
	SYMBOLS USED WITH DEMEDITEC ASSAYS	8

1 INTENDED USE

Anti-LC-1 is used for the quantitative determination of IgG antibodies to formiminotransferase-cyclodeaminase in human serum or plasma for the diagnosis of autoimmune hepatitis (AIH).

The group of primary autoimmune liver disease (PAL) comprises AIH, primary biliary cirrhosis (PBC) and primary sclerosing cholangitis (PSC). The clinical picture of PAL is in most cases not different from other chronic liver diseases. About 15% of all cases with chronic liver diseases show an autoimmune pathogenesis. Therefore, after exclusion of infectious etiology especially by viruses, the determination of different autoantibodies is recommended. Patients suffering from AIH show a variety of autoantibodies. Due to the appearance of different antibody specificities classification of AIH into different subgroups is discussed. Type I is characterized by the occurrence of antinuclear antibodies (ANA) and antibodies to smooth muscles (ASMA). ASMA recognize antigenic structure formed by polymeric f-Actin. For type II a high prevalence of antibodies to liver and kidney microsomal antigens (LKM) has been described. LKM1 antibodies recognize epitopes of cytochrome P450 IID6, a 50 kDa cytoplasmic protein found in hepatocytes and proximal tubular kidney cells. LC-1 antibodies are specific for type II hepatitis, too. The respective antigen is formiminotransferase-cyclodeaminase (FTCD) a 58 kDa metabolic enzyme located in the cytosol of liver cells. FTCD is a bifunctional protein that comprises distinct globular domains of formiminotransferase and cyclodeaminase combined by a short linker. Patients with type III autoimmune hepatitis exhibit antibodies to the soluble liver antigen (SLA). This type III is not yet fully accepted as independent subgroup for autoimmune hepatitis by the "International Hepatitis Group". In 10% of AIH cases LC-1 autoantibodies are the only liver-related circulating autoantibodies what underlines the importance of this serological parameter for the differential diagnosis of AIH.

Czaja AL. Natural history, clinical features and treatment of autoimmune hepatitis. *Semin Liver Dis* (1987) 4: 1 – 12

Johnson PJ, McFarlane IA. Meeting Report: International Autoimmune Hepatitis Group. *Hepatology* (1993) 18: 998 – 1005

Lapierre P, Hajoui O, Homberg JC, Alvarez F: Formiminotransferase cyclodeaminase is an organ-specific autoantigen recognized by sera of patients with autoimmune hepatitis. *Gastroenterol* (1999) 116: 643-649

2 PRINCIPLE OF THE TEST

Anti-LC-1 is an enzyme immunoassay for the quantitative determination of IgG antibodies to formiminotransferase-cyclodeaminase. The antibodies of the standards, controls, and diluted patient samples react with human recombinant formiminotransferase-cyclodeaminase immobilized on the solid phase of microtiter plates. The use of highly purified human recombinant formiminotransferase-cyclodeaminase guarantees the specific binding of LC-1 antibodies of the specimen under investigation. Following an incubation period of 30 min at room temperature (RT), unbound serum components are removed by a wash step. The bound IgG antibodies react specifically with anti-human-IgG conjugated to horseradish peroxidase (HRP) within the incubation period of 30 min at RT. Excessive conjugate is separated from the solid-phase immune complexes by the following wash step.

HRP converts the colorless substrate solution of 3,3',5,5'-tetramethylbenzidine (TMB) added into a blue product. The enzyme reaction is stopped by dispensing an acidic solution (HCl) into the wells after 30 min at RT turning the solution from blue to yellow. The optical density (OD) of the solution at 450 nm is directly proportional to the amount of specific antibodies bound. The standard curve is established by plotting the antibody concentrations of the standards (x-axis) and their corresponding OD values (y-axis) measured. The concentration of antibodies of the specimen is directly read off the standard curve.

3 PATIENT SAMPLES

3.1 Specimen collection and storage

Blood is taken by venipuncture. Serum is separated after clotting by centrifugation. Plasma can be used, too. Lipaemic, hemolytic or contaminated samples should not be run. Repeated freezing and thawing should be avoided. If samples are to be used for several assays, initially aliquot samples and keep at – 20 °C.

3.2 Preparation before use

Allow samples to reach room temperature prior to assay. Take care to agitate serum samples gently in order to ensure homogeneity.

Note: *Patient samples have to be diluted 1 + 100 (v/v), e.g. 10 µl sample + 1.0 ml sample diluent (C), prior to assay.*

The samples may be kept at 2 - 8 °C for up to three days. Long-term storage requires - 20 °C.

4 TEST COMPONENTS FOR 96 DETERMINATIONS

A [SORB MT]	Microtiter plate , 12 breakable strips per 8 wells coated with human recombinant formiminotransferase-cyclodeaminase	1 vacuum sealed with desiccant
B [BUF WASH 50X]	Concentrated wash buffer, 50X sufficient for 1000 ml solution	20 ml concentrate capped white
C [SAM DIL 5X]	Concentrated sample diluent, 5X	20 ml concentrate capped white
D [CONJ]	Conjugate containing anti-human-IgG- (sheep) coupled with HRP	15 ml ready for use capped blue
E [SUB TMB]	Substrate 3,3',5,5'-tetramethylbenzidine in citrate buffer containing hydrogen peroxide	15 ml ready for use capped black
F [STOP SOLN]	Stop solution 1.0 M hydrochloric acid	15 ml ready for use capped white
1 – 6 [CAL]	Standards (diluted serum) conc.: 0, 3, 10, 30, 100, 300 U/ml	1.5 ml each ready for use
P [CONTROL +]	Positive control (diluted serum) conc.: see QC data sheet	1.5 ml ready for use capped red
N [CONTROL -]	Negative control (diluted serum) conc.: see QC data sheet	1.5 ml ready for use capped green

4.1 Materials required

- micropipette 100 - 1000 µl
- micropipette 10 - 100 µl
- 8-channel wash comb with vacuum pump and waste bottle or microplate washer
- microplate reader with optical filters for 450 nm and 620 nm or 690 nm
- distilled or de-ionized water

4.2 Size and storage

Anti-LC-1 has been designed for 96 determinations. The expiry date of each component is reported on its respective label that of the complete kit on the box labels. Upon receipt, all components of the Anti-LC-1 have to be kept at 2 - 8 °C, preferably in the original kit box. After opening all kit components are stable for at least 2 months, provided proper storage.

4.3 Preparation before use

Allow all components to reach room temperature prior to use in the assay. The microtiter plate is vacuum-sealed in a foil with desiccant. The plate consists of a frame and strips with breakable wells. Allow the sealed microplate to reach room temperature before opening. Unused wells should be stored refrigerated and protected from moisture in the original cover carefully resealed. Prepare a sufficient amount of **wash solution** by diluting the concentrated wash buffer 50 times with de-ionized or distilled water. For example, dilute 1 ml of the concentrate with 49 ml of distilled water per strip.

The wash solution prepared is stable at 2 - 8 °C up to 30 days. Prepare a sufficient amount of **sample diluent** by diluting the concentrated diluent 5 times with de-ionized or distilled water. For example, dilute 10 ml of the concentrate with 40 ml of distilled water. The sample diluent prepared is stable at 2- 8 °C up to 30 days. Make sure the soak time of the wash buffer in the wells is at least 5 seconds per wash cycle. Avoid exposure of the TMB substrate solution to light!

5 ASSAY PROCEDURE

- **Dilute patient sera with sample diluent (C) 1 + 100 (v/v), e.g. 10 µl serum + 1.0 ml sample diluent (C).**
- **Avoid any time shift during pipetting of reagents and samples.**
- Bring all reagents to room temperature (18-25°C) before use. Mix gently without causing foam.
- Dispense
 - 100 µl** standards (1 - 6)
 - 100 µl** positive (P) and negative (N) control
 - 100 µl** diluted patient samples
 into the respective wells.
- Incubate **30 min** at room temperature (18-25°C).
- Decant, then wash each well **five** times using **300 µl** wash solution (made of B).
- Add **100 µl** of conjugate (D) to each well.
- Incubate **30 min** at room temperature (18-25°C).
- Decant, then wash each well **five** times using **300 µl** wash solution (made of B).
- Add **100 µl** of substrate (E) to each well.
- Incubate **30 min protected from light** at room temperature (18-25°C).
- Add **100 µl** of stop solution (F) to each well and mix gently.
- Read the OD at **450 nm** versus 620 or 690 nm within **30 min** after adding the stop solution.

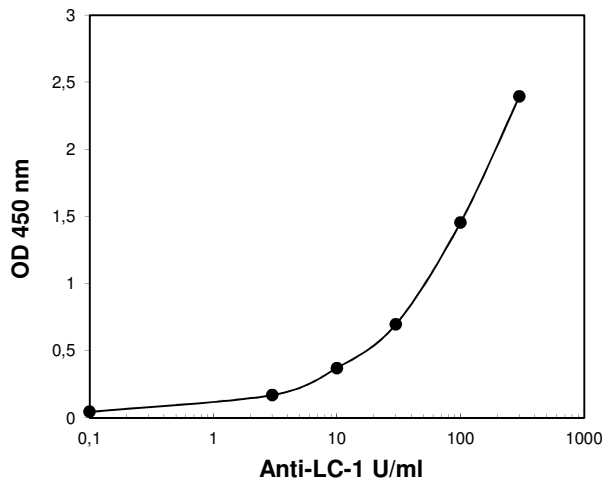
6 DATA PROCESSING

We recommend log / lin processing for best results.

The standard curve is established by plotting the mean OD-values of the standards 1 - 6 on the ordinate, y-axis, (lin. scale) versus their respective Anti-LC-1 IgG concentrations on the abscissa, x-axis, (log. scale). LC-1 antibody concentrations of the unknown samples are directly read off in U/ml against the respective OD values. Anti-LC-1 may be used also with Computer Assisted Analysis using software able to plot log/lin curves with four-parameter fit. Using the recommended dilution of 1 + 100 (v/v) for patient's sera, no correction factor is necessary, as all other components of the kit are supplied accordingly.

6.1 Example of Typical Assay Results

Well	OD (a)	OD (b)	OD (mean)	U/ml
Standard 1	0.043	0.045	0.044	0
Standard 2	0.167	0.177	0.172	3
Standard 3	0.363	0.382	0.372	10
Standard 4	0.691	0.705	0.698	30
Standard 5	1.434	1.458	1.446	100
Standard 6	2.383	2.409	2.396	300
Patient 1	0.658	0.644	0.651	26

TYPICAL STANDARD CURVE

Specimens with an OD > standard 6, should be diluted with LC-1 antibody negative serum and tested again. Results are multiplied with the dilution factor chosen. Do not use this example for interpreting results.

6.2 Test validity

The test run is valid if:

- the mean OD of the standard 1 is ≤ 0.15
- the mean OD of the standard 6 is ≥ 1.3

If the above mentioned quality criteria are not met, repeat the test and make sure that the test procedure is followed correctly (incubation times and temperatures, sample and wash buffer dilution, wash steps etc.). In case of repeated failure of the quality criteria contact your supplier.

7 REFERENCE VALUES

LC-1 antibodies	U/ml
positive	≥ 15
negative	< 15

It is recommended that each laboratory establishes its own normal and pathological reference ranges for serum Anti-LC-1 levels, as usually done for other diagnostic parameters, too. Therefore, the above mentioned reference values provide a guide only to values which might be expected.

7.1 Limitations of Method

Healthy individuals should be tested negative by the Anti-LC-1. However, LC-1 autoantibody positive apparently healthy persons do occur. Any clinical diagnosis should not be based on the results of in vitro diagnostic methods alone. Physicians are supposed to consider all clinical and laboratory findings possible to state a diagnosis.

8 CHARACTERISTIC ASSAY DATA**8.1 Calibration**

No international reference material for this parameter is available thus the assay is calibrated in arbitrary units.

8.2 Linearity

Selected positive serum samples have been tested by this assay and found to dilute linearly. However, due to the heterogeneous nature of human autoantibodies there might be sera that do not follow this rule.

8.3 Sensitivity

The analytical sensitivity of this assay was determined at 1.0 U/ml. The diagnostic sensitivity for LKM positive AIH is 48%

8.4 Specificity

No cross reactivity to other autoantigens have been found.

8.5 Precision**Intraassay variability**

sample	mean IU/ml	standard deviation	CV (%)
1	143.0	4.7	3.2
2	83.0	3.4	4.1
3	19.0	0.72	3.8

Interassay variability

sample	mean IU/ml	standard deviation	CV (%)
1	141.0	3.9	2.8
2	85.0	3.6	4.2
3	23.0	1.17	5.1











9 INCUBATION SCHEME

Dilute patients sample 10 µl serum + 1.0 ml sample diluent (C)				
1	Bring all ready for use reagents to room temperature (18-25°C) before use.			
2	Pipette	Standards (1 - 6) Control (P + N) prediluted 1 + 100 patient sera	100 µl	100 µl 100 µl
3	Incubate	30 minutes at room temperature (18-25°C)		
4	Wash	Decant, 5 x 300 µl (made of B)		
5	Pipette conjugate (D)	100 µl	100 µl	100 µl
6	Incubate	30 minutes at room temperature (18-25°C)		
7	Wash	Decant, 5 x 300 µl (made of B)		
8	Pipette substrate (E)	100 µl	100 µl	100 µl
9	Incubate protected from light	30 minutes at room temperature (18-25°C)		
10	Pipette stop solution (F)	100 µl	100 µl	100 µl
11	Measure 450 nm versus 620 (690) nm within 30 min.			

10 SAFETY PRECAUTIONS

- **This kit is for in vitro use only.** Follow the working instructions carefully. DEMEDITEC and its authorized distributors shall not be liable for damages indirectly or consequentially brought about by changing or modifying the procedure indicated. The kit should be performed by trained technical staff only.
- The expiration dates stated on the respective labels are to be observed. The same relates to the stability stated for reconstituted reagents.
- Do not use or mix reagents from different lots.
- Do not use reagents from other manufacturers.
- Avoid time shift during pipetting of reagents.
- All reagents should be kept at 2 - 8 °C before use in the original shipping container.
- Some of the reagents contain small amounts of Thimerosal (< 0.1 % w/v) and Kathon (1.0 % v/v) as preservative. They must not be swallowed or allowed to come into contact with skin or mucosa.
- Source materials derived from human body fluids or organs used in the preparation of this kit were tested and found negative for HBsAg and HIV as well as for HCV antibodies. However, no known test guarantees the absence of such viral agents. Therefore, handle all components and all patient samples as if potentially hazardous.
- Since the kit contains potentially hazardous materials, the following precautions should be observed:
 - Do not smoke, eat or drink while handling kit material,
 - Always use protective gloves,
 - Never pipette material by mouth,
 - Wipe up spills promptly, washing the affected surface thoroughly with a decontaminant.

SYMBOLS USED WITH DEMEDITEC ASSAYS

Symbol	English	Deutsch	Français	Español	Italiano
	Consult instructions for use	Gebrauchsanweisung beachten	Consulter les instructions d'utilisation	Consulte las instrucciones de uso	Consultare le istruzioni per l'uso
	European Conformity	CE-Konformitätskennzeichnung	Conformité aux normes européennes	Conformidad europea	Conformità europea
	In vitro diagnostic device	In-vitro-Diagnostikum	Usage Diagnostic in vitro	Para uso Diagnóstico in vitro	Per uso Diagnostica in vitro
	For research use only	Nur für Forschungszwecke	Seulement dans le cadre de recherches	Sólo para uso en investigación	Solo a scopo di ricerca
	Catalogue number	Katalog-Nr.	Numéro de catalogue	Número de catálogo	Numero di Catalogo
	Lot. No. / Batch code	Chargen-Nr.	Numéro de lot	Número de lote	Numero di lotto
	Contains sufficient for <n> tests/	Ausreichend für "n" Ansätze	Contenu suffisant pour "n" tests	Contenido suficiente para <n> ensayos	Contenuto sufficiente per "n" saggi
	Storage Temperature	Lagerungstemperatur	Température de conservation	Temperatura de conservación	Temperatura di conservazione
	Expiration Date	Mindesthaltbarkeitsdatum	Date limite d'utilisation	Fecha de caducidad	Data di scadenza
	Legal Manufacturer	Hersteller	Fabricant	Fabricante	Fabbricante
Distributed by	Distributor	Vertreiber	Distributeur	Distribuidor	Distributore
Content	Content	Inhalt	Conditionnement	Contenido	Contenuto
Volume/No.	Volume / No.	Volumen/Anzahl	Volume/Quantité	Volumen/Número	Volume/Quantità